


TIMING

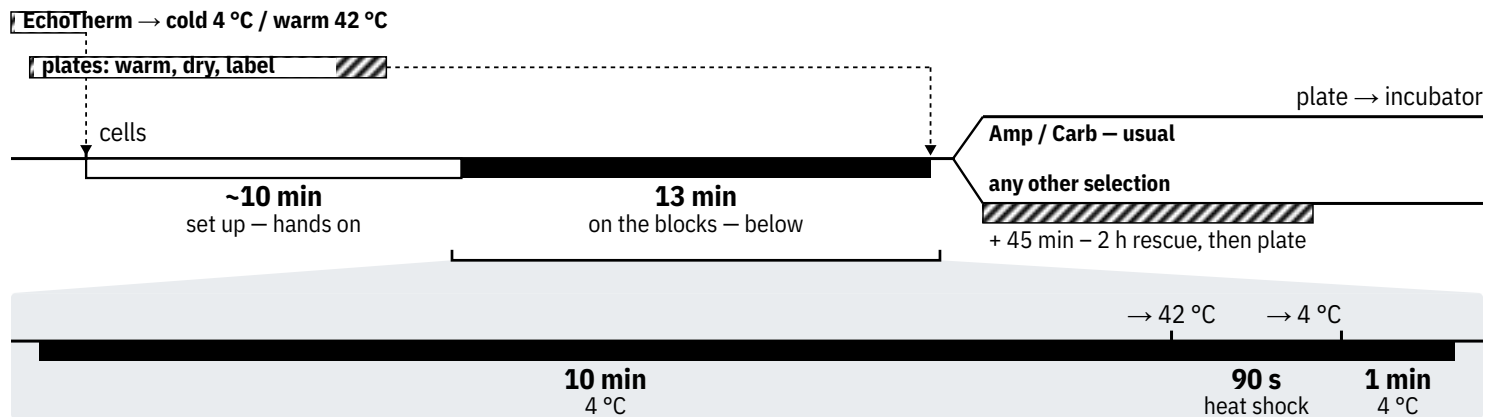

PROTOCOL

- Plates in the incubator first** — they are the slow thing. Check they carry the antibiotic(s) on your labsheet. Cold plates carry condensation and **you cannot write on a wet plate**.
- Turn on the **EchoTherm: cold block 4 °C, warm block 42 °C**. Leave both there throughout.
- Once warm and dry, **label plate bottoms**: date, initials, strain, plasmid, selection.
- Put a **100 µL** cell aliquot and the DNA tube on the **cold block**. Thaw ~30 s, add **25 µL KCM** to the cells, pipette gently. One aliquot does **3 reactions**.
- Add **40 µL** of the cell/KCM mix to the DNA. Mix gently.
 - Assumes DNA is ~20% of the total; **too much dilutes the salts**.
 - Retransforming a miniprep: **0.5 µL** into 10 µL cells.
- Run the block sequence on the timeline above**.
- Rescue only if not Amp/Carb**: 200 µL 2YT, 1.5 mL tube, shake 37 °C **45 min – 2 h. Not less, not more**.
- Plate everything, spreading with **beads**. Incubate **inverted** at 37 °C overnight, and **cancel the temperature programs**.