


KCM Heat-Shock Transformation

transformation

scan for the
full protocol



TIMING

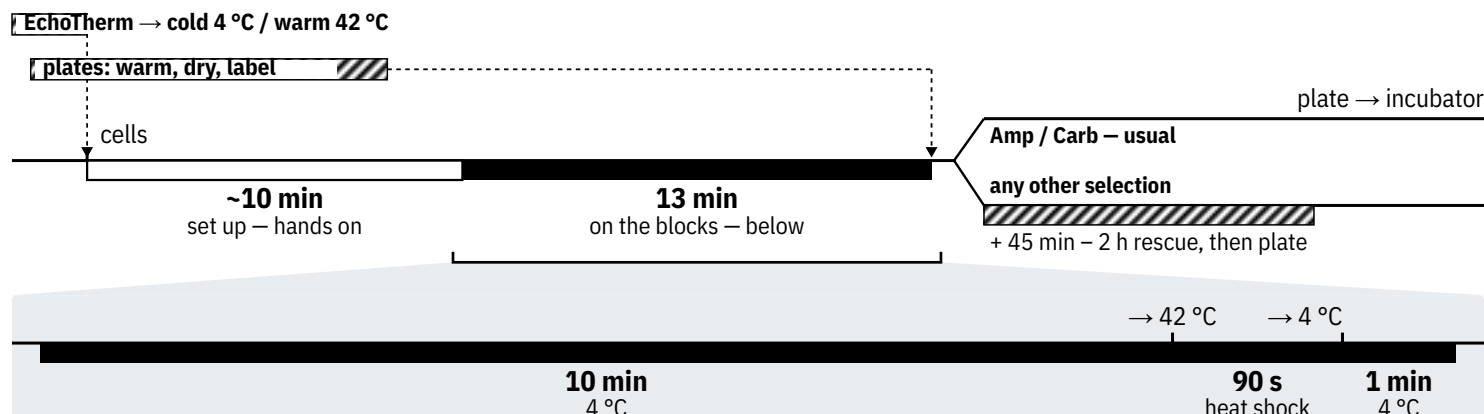

PROTOCOL

1. **Plates in the incubator first** — they are the slow thing. Check they carry the antibiotic(s) on your labsheet. Cold plates carry condensation and **you cannot write on a wet plate.**
2. Turn on the **EchoTherm: cold block 4 °C, warm block 42 °C.** Leave both there throughout.
3. Once warm and dry, **label plate bottoms:** date, initials, strain, plasmid, selection.
4. Put a **100 µL** cell aliquot and the DNA tube on the **cold block.** Thaw ~30 s, add **25 µL KCM** to the cells, pipette gently. One aliquot does **3 reactions.**
5. Add **40 µL** of the cell/KCM mix to the DNA. Mix gently.
– Assumes DNA is ~20% of the total; **too much dilutes the salts.**
Retransforming a miniprep: **0.5 µL** into 10 µL cells.
6. **Run the block sequence on the timeline above.**
7. **Rescue only if not Amp/Carb:** 200 µL 2YT, 1.5 mL tube, shake 37 °C **45 min — 2 h. Not less, not more.**
8. Plate everything, spreading with **beads.** Incubate **inverted** at 37 °C overnight, and **cancel the temperature programs.**