


KCM Heat-Shock Transformation

transformation

scan for the
full protocol



TIMING

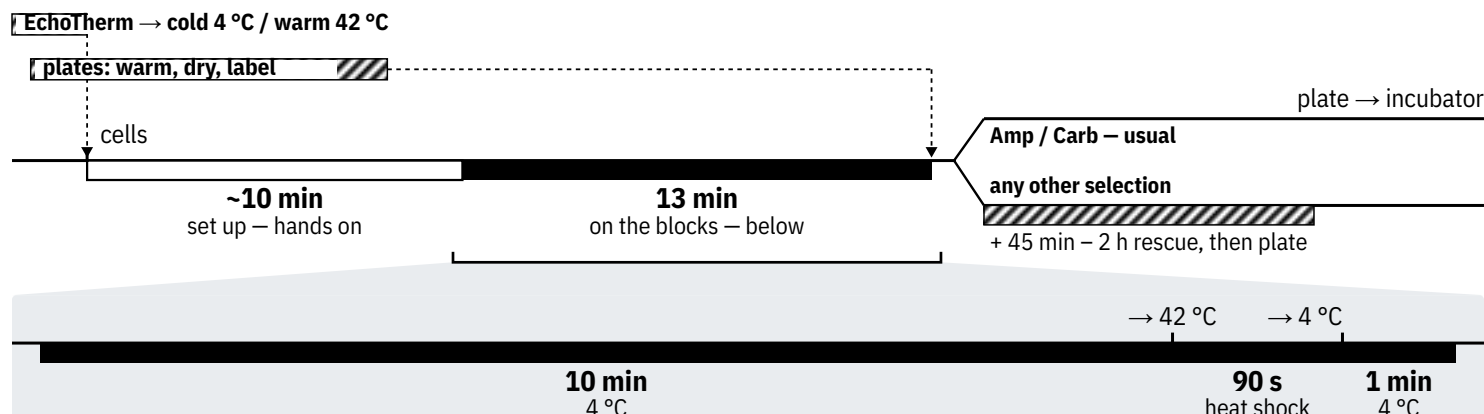

PROTOCOL

- 1. Plates in the incubator first** – they are the slow thing. Check they carry the antibiotic(s) on your labsheet. Cold plates carry condensation and **you cannot write on a wet plate.**
- Turn on the **EchoTherm: cold block 4 °C, warm block 42 °C.** Leave both there throughout.
- Once warm and dry, **label plate bottoms:** date, initials, strain, plasmid, selection.
- Put a **100 µL** cell aliquot and the DNA tube on the **cold block.** Thaw ~30 s, add **25 µL KCM** to the cells, pipette gently. One aliquot does **3 reactions.**
- Add **40 µL** of the cell/KCM mix to the DNA. Mix gently.
 - Assumes DNA is ~20% of the total; **too much dilutes the salts.**
 - Retransforming a miniprep: **0.5 µL** into 10 µL cells.
- Run the block sequence on the timeline above.**
- Rescue only if not Amp/Carb:** 200 µL 2YT, 1.5 mL tube, shake 37 °C **45 min – 2 h. Not less, not more.**
- Plate everything, spreading with **beads.** Incubate **inverted** at 37 °C overnight, and **cancel the temperature programs.**